




KEYSPLINT HARD® CLEAR

Elegoo Mars 3 Pro 4k

PRINTER SETTINGS

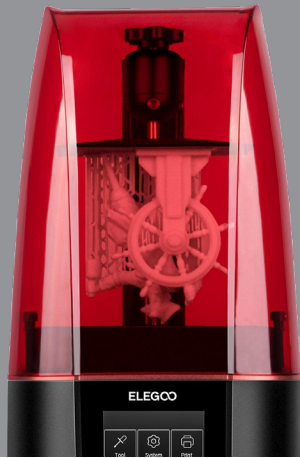

KeySplint Hard® Clear

Machine Size Scaling	X: 141.996 Y:88.709
Layer Height	0.1
Bottom Layer Count	5
Transition Layer Count	10
Exposure Time	13
Bottom Exposure Time	35
Light Off Delay*	3
Bottom Light Off Delay*	3
Bottom Lift Distance	8
Lifting Distance	8
Bottom Lift Speed	100
Lifting Speed	100
Retract Speed	150
Bottom Retract Speed	150