




KEYSPLINT SOFT®

Elegoo Mars 3 Pro 4k

PRINTER SETTINGS

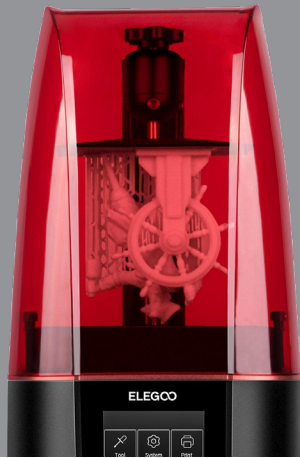

KeySplint Soft®

Machine Size Scaling	Default
Layer Height	0.1
Bottom Layer Count	5
Transition Layer Count	10
Exposure Time	13
Bottom Exposure Time	18
Light Off Delay*	3
Bottom Light Off Delay*	3
Bottom Lift Distance	8
Lifting Distance	8
Bottom Lift Speed	80
Lifting Speed	80
Retract Speed	150
Bottom Retract Speed	150