


keyprint[®]
by keystone industries

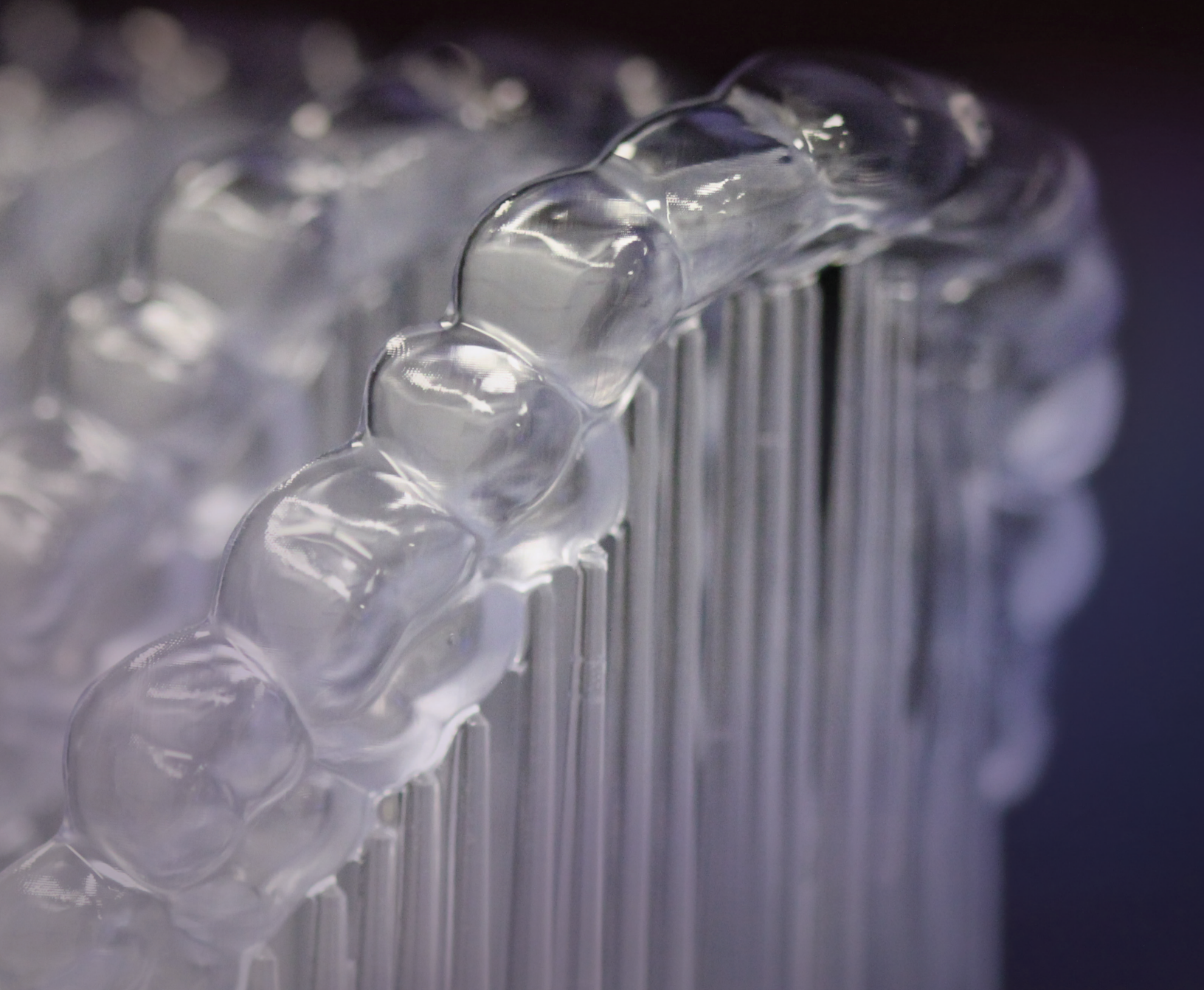

What Dental Professionals Think of KeyPrint®

“Keystone’s are the highest quality materials we have seen, and we are eagerly integrating them into our digital workflow.”
— JAY COLLINS, OWNER OF CORNERSTONE DENTAL LABS

“We’ve been receiving great feedback from dentists that KeySplint Hard® is surviving even the most severe bruxers.”
— CHRIS KABOT, VICE PRESIDENT, RESEARCH AND DEVELOPMENT AT AFFORDABLE CARE

“KeyGuide® resin is a perfect fit for our laboratory. The material prints quickly and accurately in our 3D printer and polishes to a beautiful high luster shine. Our doctors are loving it!”
— KEN KINCAID, CDT, RG, TE

“KeySplint Soft® has the unique ability to change flexural properties over the temperature range of the indicated uses. These dynamic properties add a unique characteristic that should extend the life of the splint device by preventing the degradation due to brittle fracture.”
— MARK LATTA, DMD, MS


Index	
What Dental Professionals Think of KeyPrint®	2
Our KeyPrint® Journey	4
Biocompatible & Medical Device Resins	6
KeyGuard™	8
KeySplint Soft®	10
KeySplint Hard®	14
KeyDenture Try-In™	16
KeyGuide®	18
KeyTray™	20
KeyOrtho IBT™	22
KeyPrint Workflow Validation & Open Source Philosophy	24
Laboratory Resins	28
KeyModel Ultra™	30
KeyModel®	32
KeyOrthoModel®	34
KeyMask®	36
KeyCast®	38
Accessories	40
Finishing and Polishing Materials	42
KeyPrint® Kit Accessories & Compatible Products	44
KeyPrint® Compatible	46


Our KeyPrint Journey

In an industry that is increasingly defined by the ability to innovate, adapt, and produce consistent results, Keystone Industries introduced KeyPrint in 2018, a line of precision 3D resins specifically formulated for the dental industry. Keystone has innovatively combined its unique blend of dental and photopolymer expertise to deliver the highest quality 3D printing resins to the dental market.

Entirely designed and produced in the USA, the full KeyPrint line is available throughout North America, Europe, and other global territories. KeyPrint resins are compatible with DLP, LED and LCD printers (wavelengths of 385-405 nm). Keystone is profiling 3D printers on an ongoing basis to ensure our KeyPrint resins work reliably in multiple open-source printers.

Prioritizing clinical results and patient safety, KeyPrint resins are subjected to a battery of biocompatibility tests in accordance with global regulatory standards. Each biocompatible print complies with international standards to ensure efficacy and effectiveness for the patient. KeyPrint’s robust formulations guarantee a 3-year shelf life, along with color stability.

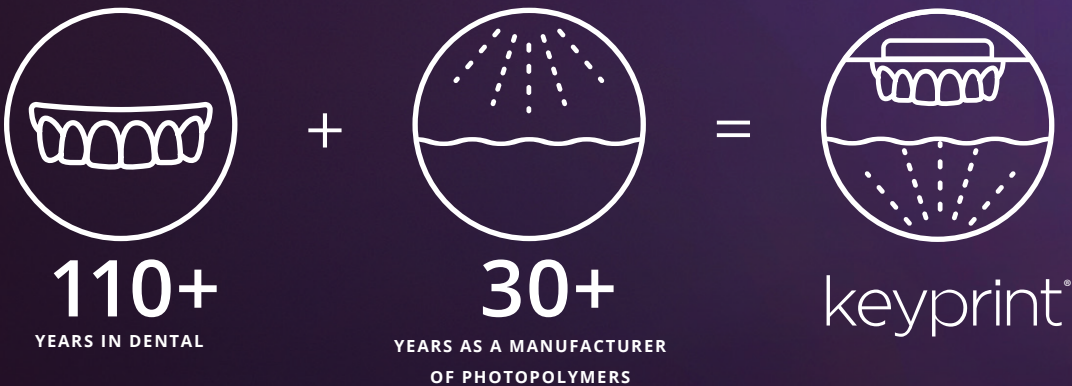

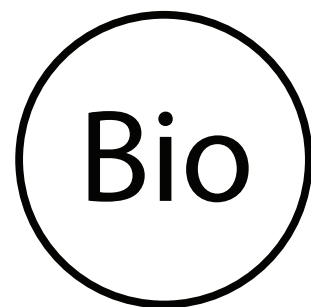

SAFE AND EFFECTIVE

Biocompatible & Medical Device Resins



WHITE

#4220200 (1kg bottle)

BLACK

#4220201 (1kg bottle)

PURPLE

#4220202 (1kg bottle)

GREEN

#4220203 (1kg bottle)


RED

#4220204 (1kg bottle)

AVAILABLE IN THE USA, CANADA, EU,
AND OTHER GLOBAL TERRITORIES




VIEW KEYGUARD
FINISHING AND
POLISHING
RECOMMENDATIONS

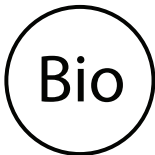


BIOCOMPATIBLE



STRONG

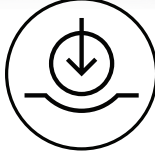

DURABLE

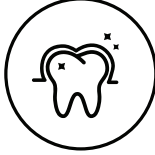

COMFORTABLE

KeyGuard™

KeyGuard™ precision dental resin produces custom-fit mouthguards for unbeatable comfort and protection!



FEATURES

- BIOCOMPATIBLE
- PRECISE FIT FOR ENHANCED COMFORT AND PERFORMANCE
- PROTECTION ONLY WHERE NEEDED WITH COMPLETE DESIGN FREEDOM
- OPTIMIZED COST EFFICIENCY
- SUPERIOR DURABILITY
- MORE SUSTAINABLE THAN TRADITIONAL METHODS

PROPERTY	METHOD	KEYGUARD
Tear Strength	ASTM D624-00	30 kN/m
Hardness	ASTM D2240-00, 25C @ 0 sec	75-90 A
Impact Strength	ISO 179-1	27.5 kJ/m2
Impact Resistance	ISO 179-1	220 J/m
Water Sorption	ASTM D570 @ 3hr	5.7%
Tensile Cyclic Hysteresis	Energy returned @ 25C Energy lost @ 25C	27% 73%
Elongation at break (strain)	ASTM D638	130%
Liquid Viscosity	Internal	900cP
Biocompatibility	ISO 10993-5, Cytotoxicity ISO 10993-10, Sensitization ISO 10993-23, Irritation	PASS

FDA 510K CLEARED | CE MARKED (CLASS IIa) | HEALTH CANADA LICENSED
AVAILABLE IN THE USA, CANADA, EU, AUS AND OTHER GLOBAL TERRITORIES

“I honestly believe this will change printed splints forever. It has unbelievable fracture resistance and mind blowing elastic memory.” — JAY PATRICK, DENTAL DEVICES BY JASON

Utilizing unique technologies, Keystone Industries has developed a one-of-a-kind 3D printing resin for the fabrication of flexible, durable dental splints and nightguards. This DLP-optimized resin is an exceptional printing material that combines the strength needed to protect teeth against disorders like bruxism, with added flexibility for better patient comfort.




VIEW OUR KEYSPLINT
SOFT VIDEO!

TINTED

- #4200005 (0.5kg bottle)
- #4220005 (1kg bottle)

CLEAR


- #4200161 (0.5kg bottle)
- #4220161 (1kg bottle)

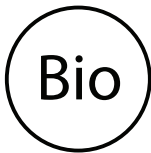

BIOCOMPATIBLE



STRONG

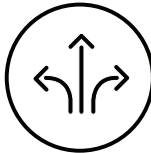

FLEXIBLE



EASY TO POLISH



EASILY CLEANED

KeySplint Soft®

Revolutionizing Night Guard and Splint Technology



FEATURES

- STRONG & FLEXIBLE; NOT BRITTLE
- BIOCOMPATIBLE MEDICAL DEVICE
- TRANSLUCENT
- UNIQUE TECHNOLOGY
- BPA, MMA & PHTHALATE FREE

PROPERTY	PROCEDURE	KEYSPLINT SOFT
Ultimate Flexural Strength	ISO 20795-2	2.6–4.4 MPa
Flexural Modulus	ASTM D790	1356
Tensile Strength	ASTM D638	52
Young’s Modulus Tensile	ASTM D638	1790
Water Sorption	ISO 20795-2	< 18 ug/mm3
Elongation at Break	ASTM D638	110%
Biocompatibility	ISO 10993-1 and USP Class VI	Pass

Doctors are loving KeySplint Soft® for their patients and their practice!

“KeySplint Soft offers unmatched durability and patient comfort. Combine that with the faster turnaround and easy replacement for night guard/splint fabrication, and you’ve got a solution that will have doctors wondering how they ever managed to practice without it.”

— DR. BRANDON OWEN, OWEN ORTHODONTICS LLC, LAKEWAY, TX

“KeySplint Soft’s unique thermoplastic features allow me to design guards that are snug but the patient loves the comfortable feel.”

— MICHAEL D SCHERER, DMD, MS, APC & MELISSA D SHOTELL, DMD, MS, SONORA, CA


CHECK OUT OUR
INDEPENDENT
WEAR STUDY!



✔ **STRONG YET FLEXIBLE; NOT BRITTLE**
Combines strength needed to protect teeth with added flexibility for better patient comfort.



✔ **EASILY CLEANED AND MAINTAINED**
Easily cleaned by the patient.

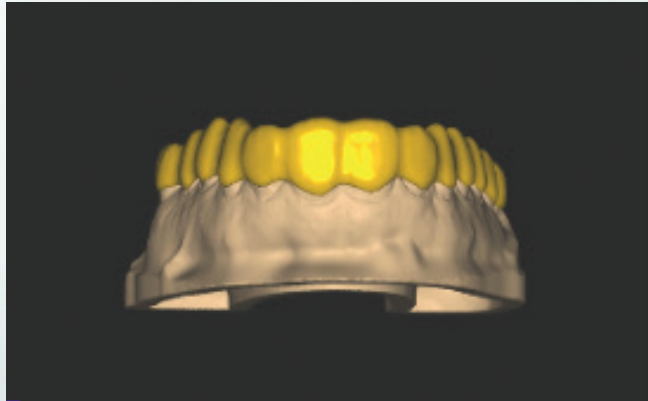

✔ **HIGHLY PRECISE DEVICES**
Precise, strong and smooth finished parts. Limited remakes and repairs can be done quickly if needed, saving time for doctors!



✔ **ELASTIC MEMORY FOR EASY PLACEMENT**
KeySplint Soft can be warmed in water for easy placement and minimal chairside adjustments, saving time for doctors!

FDA 510K CLEARED | CE MARKED (CLASS IIa) | HEALTH CANADA LICENSED
AVAILABLE IN THE USA, CANADA, EU, AUS AND OTHER GLOBAL TERRITORIES

“The KeySplint Hard material exhibited excellent wear resistance- equal to the lab processed Lucitone 199. Further the wear pattern suggests excellent wear resistance without tearing or cracking, suggesting excellent resistance to fatigue failure.” — DR. MARK A. LATTA, DMD, MS




TINTED

#4200004 (0.5kg bottle)
#4220004 (1kg bottle)

CLEAR

#4200011 (0.5kg bottle)
#4220011 (1kg bottle)

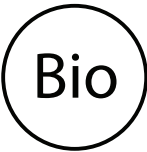


BIOCOMPATIBLE




DURABLE



EASY TO POLISH

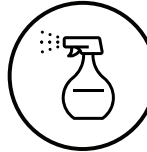

EASILY CLEANED

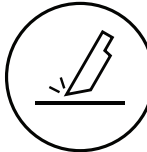

ABRASION RESISTANT

KeySplint Hard®

The newest member of the KeySplint family, KeySplint Hard is designed for rigid dental splints. KeySplint Hard provides therapy for bruxism and cases where tooth immobilization is needed.



- FEATURES
- BIOCOMPATIBLE MEDICAL DEVICE
 - UNITED STATES FDA 510(K) CLEARED
 - HEALTH CANADA LICENSED
 - CE MARKED (CLASS IIA)
 - HARD FOR RIGID SPLINTS
 - BPA, MMA & PHTHALATE FREE



PROPERTY	PROCEDURE	KEYSPLINT HARD
Flexural Strength (MPa)	ISO 20795-2	60-65 MPa
Flexural Modulus (MPa)	ISO 20795-2	1510-1600 MPa
Water Sorption	ISO 20795-2	18 µg/mm ³
Elongation at Break	ASTM D638	9%
Biocompatibility	ISO 10993-1	Pass

Class I

AVAILABLE IN THE USA, CANADA, EU, AUS
AND OTHER GLOBAL TERRITORIES

- A1

#4200401 (0.5kg bottle)
#4220401 (1kg bottle)
- B1

#4200406 (0.5kg bottle)
#4220406 (1kg bottle)
- BL4


#4200420 (0.5kg bottle)
#4220420 (1kg bottle)



BIOCOMPATIBLE



COST
EFFECTIVE

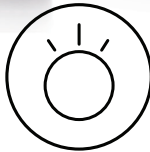

EXCEPTIONAL
SURFACE

KeyDenture Try-In™



KeyDenture Try-In takes the uncertainty out of the denture fitting process. Designed for precision and efficiency, KeyDenture Try-In is a photopolymer resin indicated for 3D printing of biocompatible, short-term trial dentures, intended for evaluation of fit, occlusion, esthetics, and phonetics during try-in appointments.



- FEATURES**
- COST EFFECTIVE FOR TRY-IN DENTURES
 - EXCEPTIONAL SURFACE QUALITY
 - EASY TO CHARACTERIZE
 - POLISHES TO HIGH LUSTER
 - BIOCOMPATIBLE MEDICAL DEVICE
 - WIDE VARIETY OF KEYPRINT COMPATIBLE WORKFLOWS
 - BPA, MMA & PHTHALATE FREE

PROPERTY	PROCEDURE	KEYDENTURE TRY-IN
Flexural Strength (MPa)	ASTM D790	56-68 MPa
Flexural Modulus (MPa)	ASTM D790	1500-1750 MPa
Water Sorption	ISO 20795-1	9-13 µg/mm3
Elongation at Break	ASTM D638	30-43%
Biocompatibility	ISO 10993-10	Pass

AVAILABLE IN THE USA, CANADA, EU, AUS AND OTHER GLOBAL TERRITORIES


“When we tested this material in the beginning, I was surprised. It prints really well and the fit was very precise. It was really easy to use. We are definitely impressed with it. We will no longer be outsourcing our surgical guides.” — JESSICA BIRRELL, CDT

Class I

- #4200001 (0.5kg bottle)
- #4220001 (1kg bottle)



VIEW OUR
STERILIZATION
PROTOCOL!

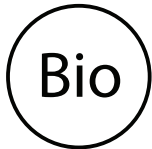


BIOCOMPATIBLE



STRONG




EASY TO POLISH

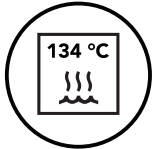

AUTOCCLAVABLE



KeyGuide®

KeyGuide is ideal for fabricating transparent and stable surgical guides, which enable the placement of implants at a precise angle and depth.

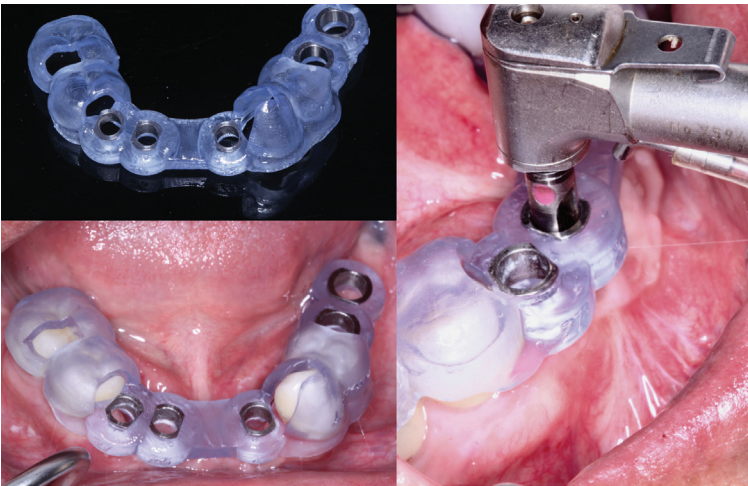

Photos courtesy of August De Oliveria, DDS.


- FEATURES
- STRONG
 - BIOCOMPATIBLE MEDICAL DEVICE
 - EASY TO POLISH
 - AUTOCLAVABLE AND COLOR STABLE
 - BPA, MMA & PHTHALATE FREE

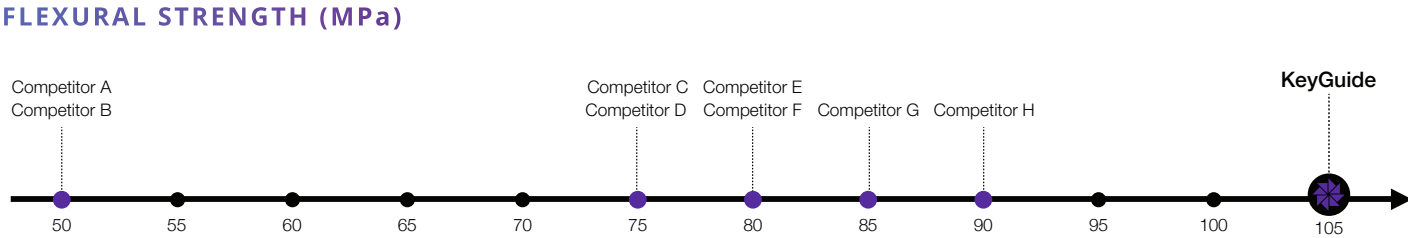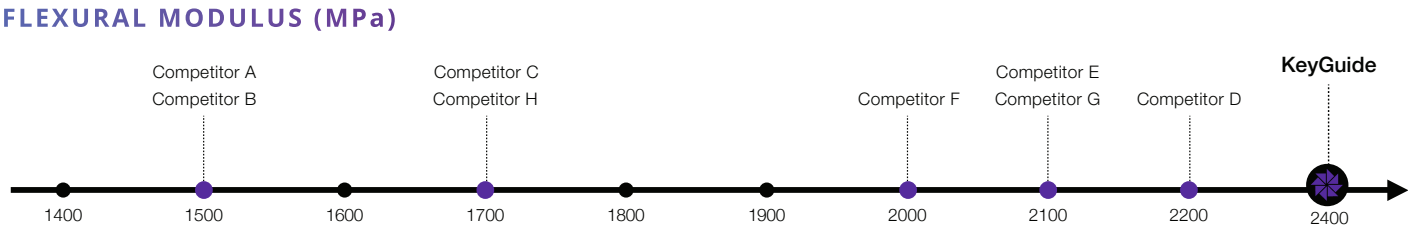


AVAILABLE IN THE USA, CANADA, EU, AUS AND OTHER GLOBAL TERRITORIES

Class I

- #4200007 (0.5kg bottle)
- #4220007 (1kg bottle)

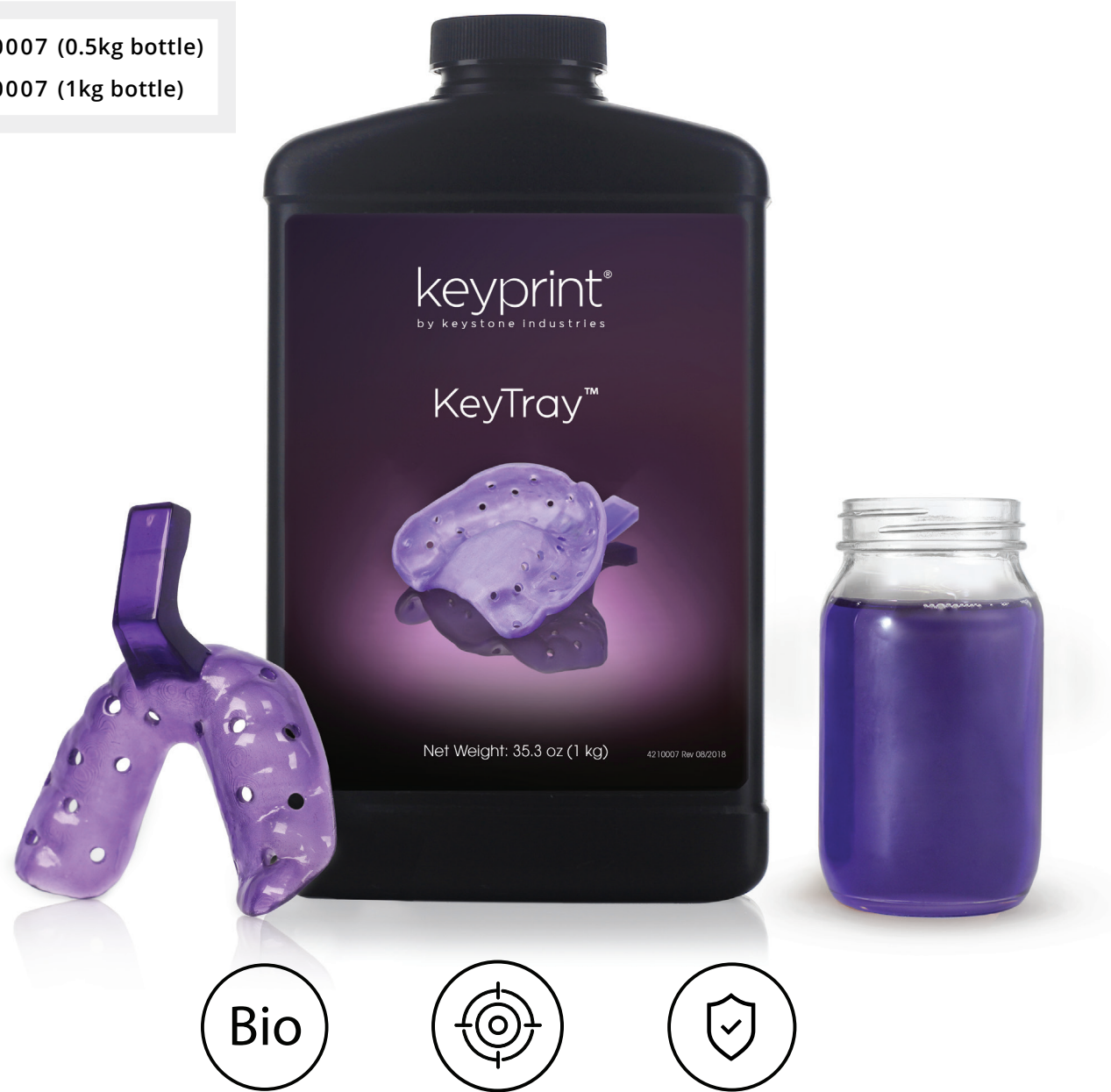

BIOCOMPATIBLE

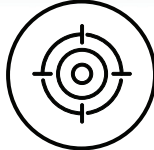

IMPRESSION
ACCURACY



STRONG

KeyTray™

KeyTray is a strong, biocompatible (Class I) 3D printing resin designed to create customized, individual impression trays quickly and with precision. The material is strong and durable to withstand the forces of taking a patient impression and removing the tray from the oral cavity. It is compatible with all types of impression material.



FEATURES

- BIOCOMPATIBLE MEDICAL DEVICE
- NO PRELIMINARY CAST REQUIRED
- SAVES TIME BY ELIMINATING STONE MODELS
- IMPROVED IMPRESSION ACCURACY
- COMPOUND WAXES AND BORDER MOLDING MATERIALS WILL ADHERE TO TRAY
- BPA, MMA & PHTHALATE FREE

PROPERTY	PROCEDURE	KEYTRAY
Flexural Strength (MPa)	ASTM D790	84.7
Flexural Modulus (MPa)	ASTM D790	1913
Tensile Strength (MPa)	ASTM D638	62.1
Tensile Young's Modulus (MPa)	ASTM D638	2056
Elongation at Break (%)	ASTMD638	26.4
Shore D Hardness		86.5 D

AVAILABLE IN THE USA, CANADA, EU, AUS AND OTHER GLOBAL TERRITORIES

“You can print these things all day long and not smell anything, the odorless feature is great. This is also the easiest to peel away of any of the different resins that we’ve ever tried for the indirect bond.” — **BRANDON OWEN DDS, MS**

Class I

- #4200008 (0.5kg bottle)
- #4220008 (1kg bottle)



VIEW OUR
KEYORTHO IBT VIDEO!



BIOCOMPATIBLE




PRECISE



FLEXIBLE



ODORLESS

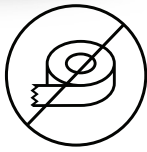

NON-STICK

KeyOrtho IBT™


KeyOrtho IBT is ideal for printing orthodontic transfer trays for positioning brackets. KeyOrtho IBT’s translucency allows for precise positioning before and during bracket placement. Its high tensile strength allows for easy application of brackets in a single step, reducing chair time to minutes. It is flexible for easy removal, and works with all general bracket systems. It is also odorless and tasteless for patient comfort.



FEATURES


- BIOCOMPATIBLE MEDICAL DEVICE
- ODORLESS & TASTELESS
- ADHESIVES RELEASE EASILY FROM TRAY
- TRANSLUCENT FOR EASY POSITIONING
- PRECISE
- FLEXIBLE
- BPA, MMA & PHTHALATE FREE

PROPERTY	PROCEDURE	KEYORTHO IBT
Tensile Strength (MPa)	ASTM D638	8.0 - 10.5
Tensile Young’s Modulus (MPa)	ASTM D638	15.5 - 31.0
Elongation at Break (%)	ASTM D638	> 130%
Biocompatibility	ISO 10993-5	Pass
Viscosity @ 25°C (cP)	ASTM D7867	<1000
Orthodontic Adhesive Release		Pass

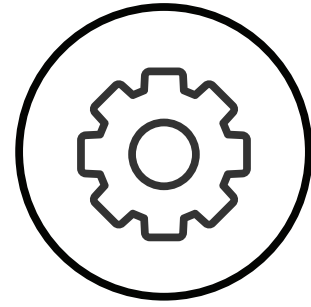

KEYPRINT WORKFLOW VALIDATION

Consumer Confidence in Safety & Efficacy

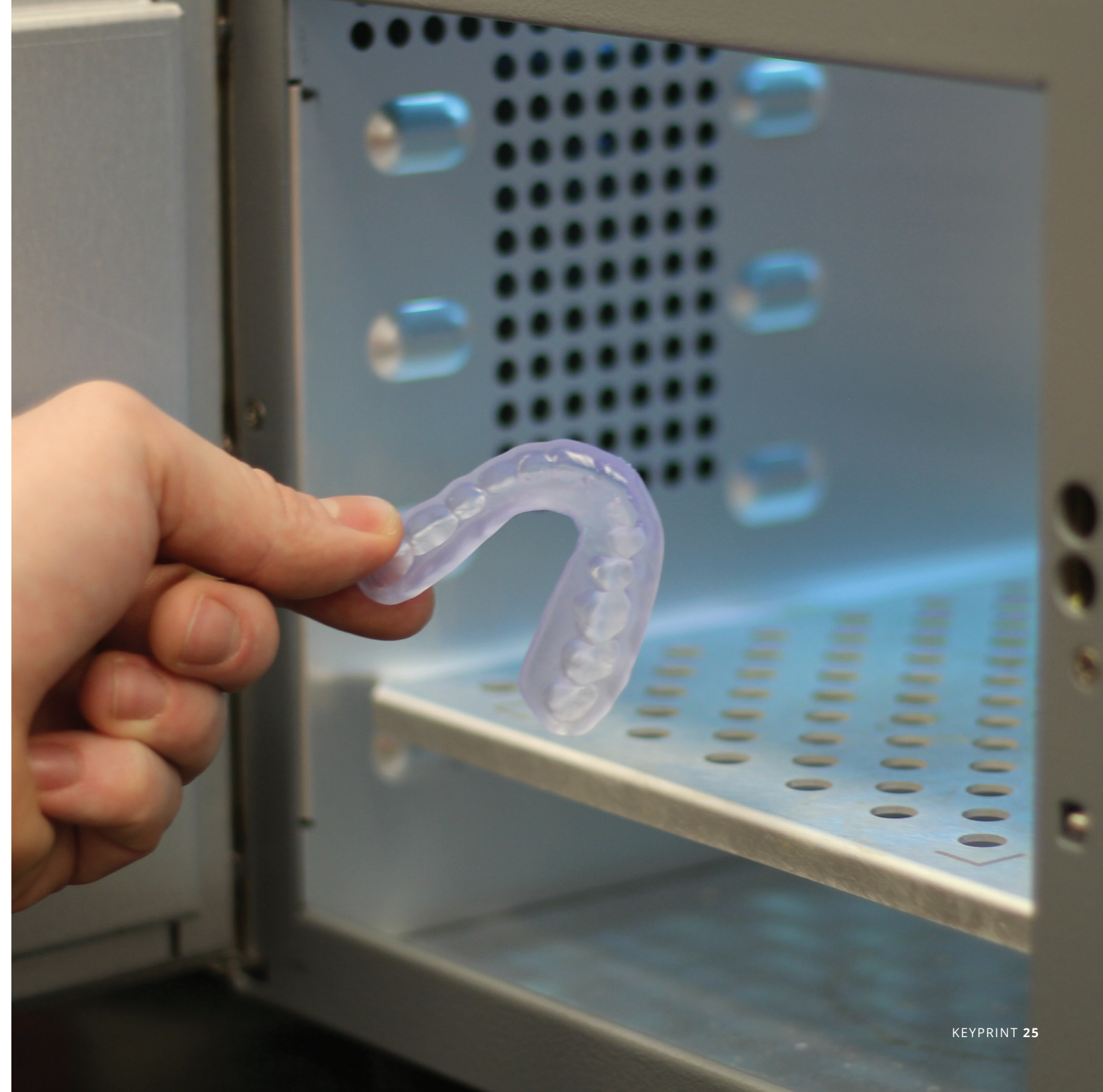

1
DESIGN



OUR OPEN SOURCE PHILOSOPHY

Keystone is always profiling new printers to ensure that KeyPrint resins work reliably in as many open-source printers as possible. We can work with you to optimize your printer settings to work with KeyPrint, but Keystone does not guarantee high-quality print results from non-validated printers.


KeyPrint resins are designed to work in 385-405 nm DLP, LED and LCD printers to best yield detailed parts faster and with superior accuracy. for customers looking to purchase resins designed for SLA and other printer technologies, Keystone is willing and able to discuss developing specialized resins for your system.

DESIGN

We work with all open source dental design software:

- ✓ 3Shape
- ✓ Dental Wings
- ✓ Exocad
- ✓ D3Tool
- ✓ Blue Sky Bio
- ✓ CADflow.ai

2
PRINT

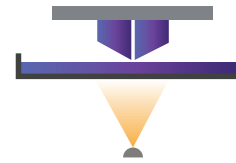

PRINTER VALIDATIONS

Choosing one of Keystone's validated partner printers guarantees:

- ✓ Printing settings tested and approved by Keystone's Technical Team
- ✓ Optimal printing accuracy for all KeyPrint Resins
- ✓ Passing biocompatibility checks on medical device resins
- ✓ Technical support from both the Printer Partner and Keystone



VIEW KEYSTONE'S PRINTER
VALIDATIONS

3
POST WASH



POST WASH

Following Keystone's validated post wash protocol is critical to a successful 3D printing digital workflow. All KeyPrint resins require post wash.

Keystone's Recommended post wash includes the following:

- ✓ Optimize™ IPA 97%+ purity, following the 3 bath protocol listed in the material's Instructions For Use (IFU)
- ✓ Submerged agitation is necessary to remove all residual resin. A soft brush is recommended to help with this process
- ✓ 3D Printed parts should be submerged for **no more than 5 minutes total** in the solvent
- ✓ It is recommended to use compressed air after post wash to ensure all residual resin and IPA is removed from the part
- ✓ Allow the KeyPrint parts to dry (covered) for at least 30 minutes before curing

4
POST CURE




POST CURE

Similar to green-state zirconia, 3D printed parts require additional treatment after printing to ensure they have :

- ✓ Biocompatibility for intraoral use
- ✓ Color stability
- ✓ Ideal physical properties

All KeyPrint resins require UV Curing after post wash. Keystone publishes the validated post cure units and cure times to ensure the indicated devices are safe and effective.

NOTE: Nitrogen is not needed

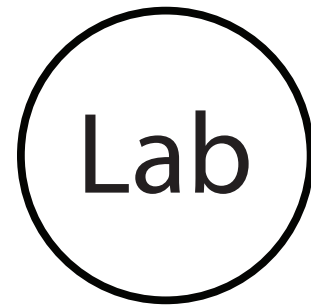

ACCURATE, RELIABLE, AND CONSISTENT

Laboratory Resins



IVORY


#4220009 (1kg bottle)
#4240009 (4kg bottle)

LIGHT GREY


#4220010 (1kg bottle)
#4240010 (4kg bottle)




VIEW OUR KEYMODEL
ULTRA VIDEO!

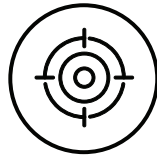

FLAWLESS
SURFACE DETAIL

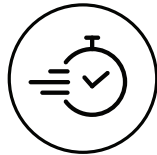

ULTRA-FAST
TECHNOLOGY

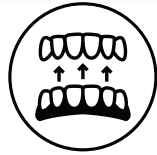

THERMOFORMING
RELEASE AGENT

KeyModel Ultra™

KeyModel Ultra is a material designed for 3D printing of dental and orthodontic models.



IDEAL FOR THERMOFORMING ALIGNERS

FEATURES


- HIGH-SPEED PRINT FORMULATION
- INTEGRATED RELEASE AGENT FOR EASY THERMOFORM SEPARATION
- TWO COLORS OPTIMIZED FOR SHARP DETAIL
- SMOOTH SURFACE FINISH
- STABILITY IN PRODUCTION WORKFLOW
- INTEGRATED WITH A WIDE RANGE OF 3D PRINTERS

PROPERTY	PROCEDURE	KEYMODEL ULTRA
Maximum Tensile Strength (MPa)	ASTM D638	50
Tensile Modulus (MPa)	ASTM D638	1700
Elongation at Break (%)	ASTM D638	51
Maximum Flexural Strength	ASTM D790	70
Flexural Modulus	ASTM D790	1940
Viscosity @ 25°C (cP)	ASTM D7867	500-600


#4200006 (0.5kg bottle)
#4220006 (1kg bottle)



ACCURATE

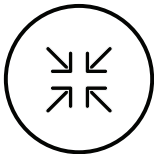

LOW SHRINKAGE

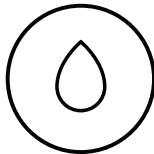

EASY VISIBILITY

KeyModel[®]

KeyModel produces highly precise and accurate dental models. Its low shrinkage formulation makes it an ideal resin for crown, bridge, and implant models, orthodontic clear aligner models, and diagnostic modeling.


- FEATURES**
- CARVABLE WITHOUT CHIPPING
 - COLORED FOR EASY VISIBILITY OF MARGINS AND TOOTH ANATOMY



OPTIMAL SCANNING
Models have smooth, hard surfaces for optimal scanning



ACCURATE & PRECISE
Highly accurate detail in printing

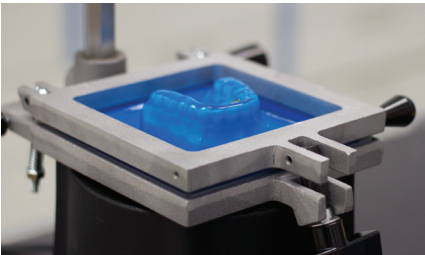

PAIRS WITH: PRO-FORM[®]
Check out our line of Pro-form thermoforming materials at dental.keystoneindustries.com

PROPERTY	PROCEDURE	KEYMODEL
Flexural Strength (MPa)	ASTM D790	79.6
Flexural Modulus (MPa)	ASTM D790	2293
Tensile Strength	ASTM D638	51
Young's Modulus (MPa)	ASTM D638	1850
Elongation at Break (%)	ASTM D638	5.60%
Shore Hardness	Scale D	86

#4200003 (0.5kg bottle)
#4220003 (1kg bottle)



KeyOrthoModel[®]


KeyOrthoModel is designed for high-speed and accurate prints. Perfect for thermoforming aligners and orthodontic devices.




- FEATURES**
- OPTIMIZED TO PRINT @ 100 AND 200µm LAYER THICKNESS FOR HIGH THROUGHPUT
 - PRINTS IN 1/2 THE TIME OF KEYMODEL AT SAME LAYER THICKNESS
 - PAIRS WITH PRO-FORM[®]
 - PERFECT FOR MAKING ALIGNERS
 - WITHSTANDS THERMOFORMING TEMPERATURES

PROPERTY	PROCEDURE	KEYORTHOMODEL
Flexural Strength (MPa)	ASTM D790	79.6
Flexural Modulus (MPa)	ASTM D790	2293
Tensile Strength	ASTM D638	51
Young's Modulus (MPa)	ASTM D638	1850
Elongation at Break (%)	ASTM D638	5.60%
Shore Hardness	Scale D	86


#4220002 (0.5kg bottle)

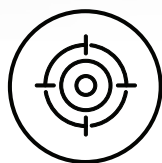

ACCURATE

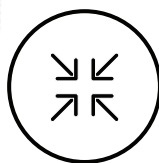

LOW SHRINKAGE

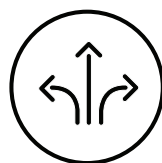

FLEXIBLE



SIMULATES
GUM TISSUE

KeyMask®

Flexible yet firm gingival mask that simulates the look and feel of gum tissue.
Incorporate KeyMask into your digital restoration planning.



FLEXIBLE YET FIRM

This gingival mask simulates the look and feel of gum tissue



EASILY TRIMMED

Creates a proper emergence profile for pontic shape


PROPERTY	PROCEDURE	KEYMASK
Tensile Strength (MPa)	ASTM D638	2
Young's Modulus (MPa)	ASTM D638	4.6
Elongation at Break (%)	ASTM D638	41%


“KeyCast prints well, with excellent detail needed for RPDs. It burnouts thoroughly and leaves no discernible residual ash. Once properly post-cured, it should slide into your existing equipment workflow for burnout and casting.” — **BY KEN KINCAID CDT,TE**

#4200000 (0.5kg bottle)
#4220000 (1kg bottle)



**ACCURATE**

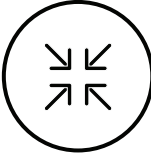**LOW SHRINKAGE**

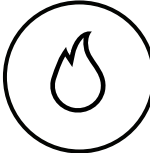
**ASH-FREE
BURNOUT**

KeyCast[®]

KeyCast’s low shrinkage formula controls thermal expansion to yield porous-free and accurate castings. KeyCast is compatible with rapid burnout and gradual temperature ramp workflows, with no residual ash. It is used with chrome, cobalt, gold, and silver.

PROPERTY	PROCEDURE	KEYCAST
Flexural Strength (MPa)	ASTM D790	15
Flexural Modulus (MPa)	ASTM D790	340
Tensile Strength (MPa)	ASTM D638	9.1
Young’s Modulus (MPa)	ASTM D638	220
Elongation at Break (%)	ASTM D638	16%
CTE 0-150 °C (32-302 °F)	ASTM E831	168 x 10-6m/(m. °C)
CTE 150-250 °C (302-482 °F)	ASTM E831	64 x 10-6m/(m. °C)
CTE >250 °C (>482 °F)	ASTM E831	<0m/(m. °C)
Shrinkage	ASTM D792	<2%
Burnout		<0.05%

PAIRS PERFECTLY WITH:

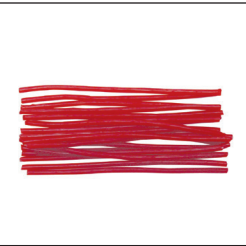

Sprue Wax




Microtorch



Green Flexiboles

Check out our full line of dental casting rings, waxes, mixing bowls and other accessories at dental.keystoneindustries.com.

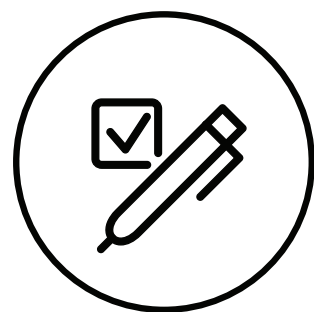

Accessories



Finishing and Polishing Materials

ITEM NAME	ITEM NUMBER	PHOTO
KeyPolish™ Kit	#1009420	
Mini Scotch Brite™	#1670092	
400 Grit Blue Disc	#1670099	
Course Grit Pumice*	#1700010	
Muslin Buffs	#1180190 for 5" X 45 #1180170 for 4" X 42	
Tripoli Polishing Compound	#1660030	
Beige Polishing Paste (100g)	#1670175	
High-Shine Cotton Buff	#1670140 or #1670510	
High-Shine Polistar Pink (150g)	#1670546	

* sold separately from kit



VIEW KEYSPLINT
FINISHING AND POLISHING
RECOMMENDATIONS



VIEW OUR
POLISHING VIDEO



BEFORE POLISHING



AFTER POLISHING

KeyPrint® Kit Accessories



DENTURE 3D BASKETS
#1220207



MICROFIBER CLOTH
#1220201



PROFESSIONAL SPRUE CUTTER
#1220200



1 1/2" STIFF PUTTY KNIFE
#1220210



1 1/2" PRO CHISEL
#1220209



WASH BOTTLE 250 ML.
#1655413



SMALL BOWL SCRAPER
#1220202



TEKTON Mallet 8 OZ.
#1220208



1 1/2" BRISTLE BRUSH
#1220205



2 1/2" BRISTLE BRUSH
#1220206



SCRATCH .016 NYLON BRUSH
#1220203



6" POINT CURVED TWEEZERS
#1220204



ALL PURPOSE UTILITY SCISSORS
#9139720

#4220010K Accessories Kit





KEYPRINT COMPATIBLE

Keystone provides high-quality materials, workflow efficiencies, and a commitment to open partnerships. We are committed to expanding our KeyPrint Compatible portfolio, providing an unmatched selection of printers and curing units to select from. Here's what it means to be KeyPrint Compatible:

- 15 UNIQUE RESINS
- 46 VALIDATED PRINTERS & COUNTING!
- 31 VALIDATED CURING UNITS
- 3,500+ PRINT WORKFLOWS VERIFIED
- 325+ PASSING BIOCOMPATIBILITY TESTS, INCLUDING 50 FOR THE KEYSPLINT FAMILY



[VIEW OUR
LATEST PRINTER
VALIDATIONS](#)

Keystone has validated the following workflows as KeyPrint Compatible with select resins. All of our resins are compatible across DLP / LCD / LED from 385-405nm.

Carbon

rapidshape

SprintRay

SHINING 3D

MICROLAY
SLA SYSTEMS

FLASHFORGE
3D PRINTER

NewPro3D

LuxCreo
Digital 3D Production

ivoclar

B9Creations

ASIGA

ACKURETTA

stratasys

HEYGEARS

MII CRAFT

UNIZ

PLANMECA

RAYSHAPE

zortrax

ANYCUBIC

nEXA3D

Desktop Health

phrozen

W2P

EPAX 3D

DEKEMA

sisma

Photo centric

ELEGOO

PRECISION 3D RESINS

HEADQUARTERS

480 South Democrat Road,
Gibbstown, NJ 08027 USA
1.800.333.3131

INTERNATIONAL

Stockholzstr. 11
D-78224 Singen, Germany
+49 77 31 91 21 01
sales@keystoneind.com

KEYPRINT.COM



JOIN OUR
KEYPRINT
COMMUNITY ON
FACEBOOK!

FDA registered and ISO 13485 certified

For more information, please call 1.800.333.3131. To view a full list of our validated printers and post-cure units, please visit: keyprint.com/compatible

REV 1 4/2024